



PathStand™ Stand Alone Grossing Solutions

Telepathology and Image Archiving for Pathology



PathStand™ | Routine

Complete stand-alone Imaging system. Supports in-process documentation.



PathStand Key Features

- 20 MP Image Capture
- 720p Live Video
- 34X zoom
- Metal Construction
- All-in-one Computer
- LED Lights
- Washable Keyboard and Mouse
- Bar Code Scanner
- Footswitch
- Floating position background board
- Imaging Glass
- Custom Color Backgrounds



Superior Image Quality


The PathStand provides HD resolution across a wide range of magnification by utilizing a highly corrected 34X zoom lens. The large zoom range accommodates large and small specimens as well as small details on large specimens.

- Max FOV: 86 cm
- Min FOV: 2.5 cm
- 720p HD Resolution

1X Zoom



8.33X Zoom

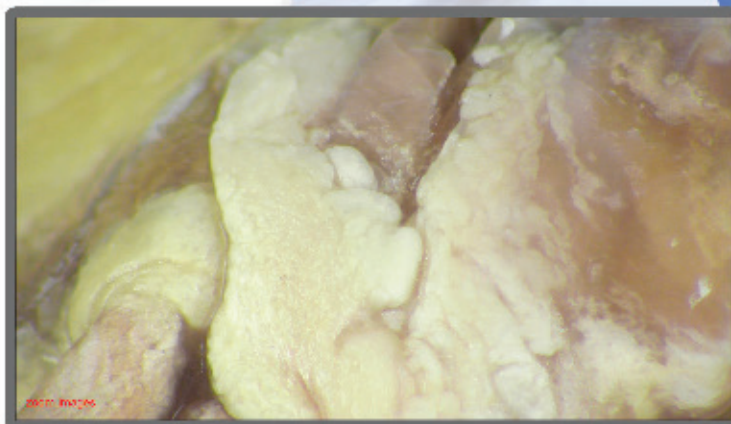

34X
ZOOM


PathStand™ | Comprehensive

Fully integrated and completely scalable - integrating cameras, audio, barcode scanners and software. Our pathology imaging suite is designed to empower you from scanning to case sign-out.



HD Video

Crisp, real time 720p HD video provides live viewing of your specimen for both local and remote telepathology users ensuring efficient and effective workflow.

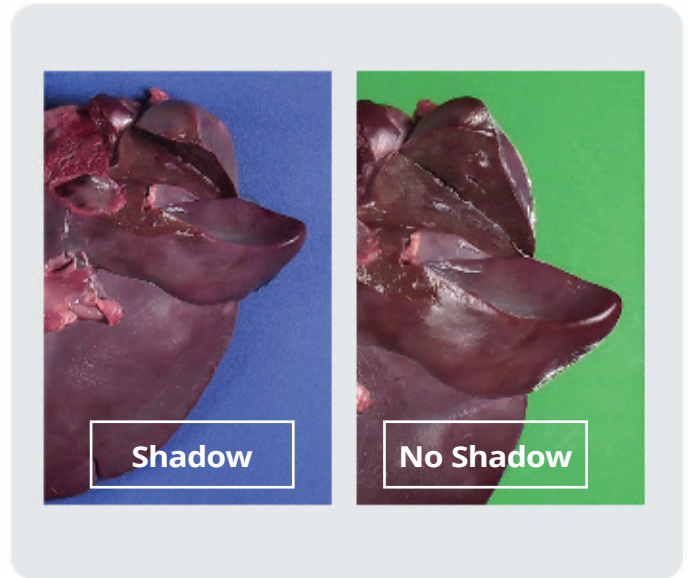

No Shadow Glass

PathStand includes both a solid board that provides a traditional imaging background and a glass table to facilitate shadowless imaging by suspending your sample above the background. Publication quality images are a snap.



>12 Year Life!

LED Lighting

Great images result from great lighting and the PathStand provides adjustable, color balanced, LED light stalks. These lights are bright, energy efficient and have a 50,000 hour life (>12 years at 70 hrs per week). Right and left banks are individually switched to provide directional lighting.



Lab Hardened


Labs are harsh environments for instruments and PathStand was designed to survive the rigors of cleaning and disinfecting agents. Components are powder coated, anodized aluminum or stainless steel to ensure carefree operation over the system life.

PathStand™ | Easy to Learn - Easy to Use

PathSuite provides the user with a simple icon driven user interface with one click access to essential controls. This makes PathSuite intuitive to learn and a pleasure to use. PathSuite supports both gross and microscopy imaging with a common user interface.

Gross Imaging

Labels in the screenshot include: Capture, Zoom, Brightness, Focus, White Balance, Full Screen, Magnifier, Teleconferencing, Setup, Case, Case Management Tools, Film Strip, Image, and Foot Pedal.

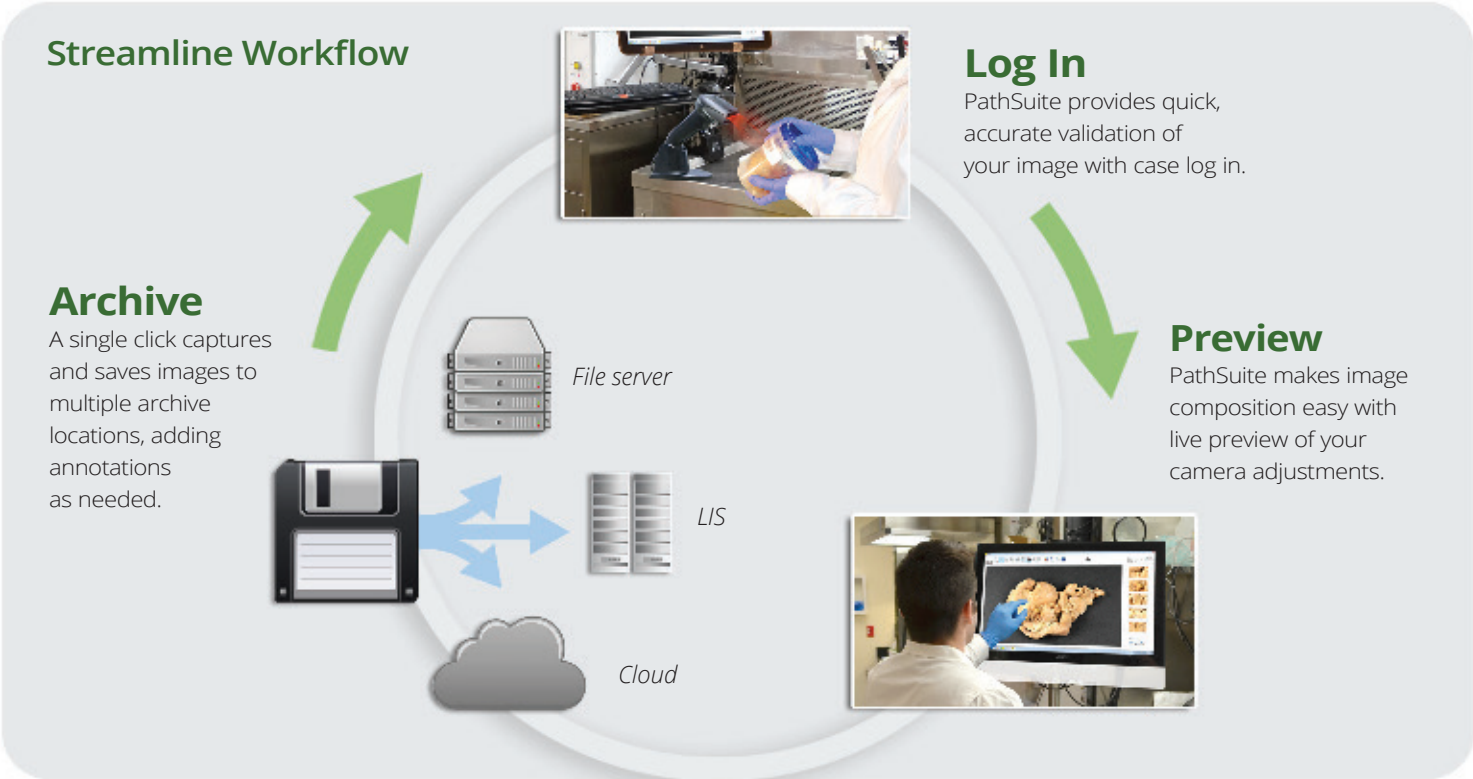

PathStand™ | Optimized Interface

Measurement and Annotate

PathSuite provides a full suite of annotation and measurement tools for detailing your images without transferring images between programs or computers.



Auto Imprints

Calibration marks and case numbers add visual reference to your image.



Measurement

PathSuite's pre-calibrated images allow for easy measurement of length, area & angle.



Custom Annotations

A full set of text and drawing tools ensure fully detailed images.



Style

Annotation style preferences are easily set or changed with control over numbering protocol, font, size, color, fill, etc...



Case Management Tools

Easily select, duplicate and move images to support your "after hours" and case management work.

Telepathology


PathSuite provides simple, secure connections and fully interactive live consults in a practical and affordable format.

Connect



Colleague meetings are initiated by selecting the team member from the connection menu. The colleague receives an invitation prompt on their screen and joins the meeting with a simple click. Non-colleague attendees are easily added to meetings via e-mail invites. Once accepted, meetings are secured using guest passwords and data encryption, supporting HIPPA requirements.

Consult

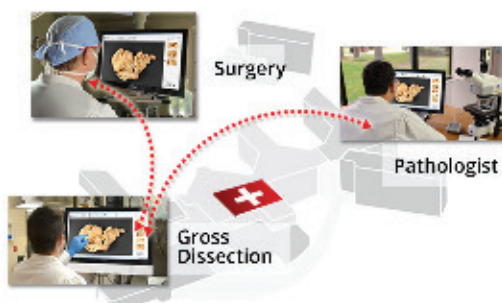


PathSuite ensures a rich collaborative consult allowing colleagues to easily clarify orientation, discuss dissection and verify margins.

Document



Discussion images are easily saved to the case with a click of a button preserving results for later reference and follow up.

Save Time

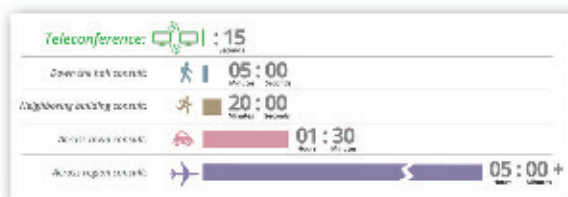

PathSuite telepathology solutions save you time on each trip you make whether it be down the hall, down the street or across the world. It eliminates the time and frustration associated with the trip and the unaccounted for impromptu interactions that occur with each trip. Telepathology lets you connect, consult and return to work without leaving your desk, it's that simple.

PathStand™ | Practical

Return on Your Investment

Our install base has garnered dramatic time savings on image acquisition across multiple pathology work tasks resulting in impressively short ROI periods, followed by ongoing savings year after year.



Routine Gross Dissection



Time Saved per Specimen **10 Min** × Specimens per Day **12** = Time Saved per Day **2 Hrs** or Time Saved per Year **520 Hrs**

Increase Consultations

Increasing intraoperative and needle biopsy adequacy consultations provides savings for the organization and the physician by eliminating unreimbursed secondary procedures and increasing patient satisfaction. Pathologists and Pathologist Assistants benefit from the increased efficiency of providing their consults virtually, avoiding the unproductive trip otherwise required.

Bundled Rate Programs: Planning for Performance-Based Reimbursement

Changing Healthcare Design



Performance-Based Bundled Payments

2012 Fall Managed Care Forum
November 5-8, 2012

Presented by **David A. Gregory, FACIC**


Professional Documentation

Consistent documentation of your work provides the foundation for quality assurance and process improvement. It is also the basis of how you are perceived by your patients, referring physicians and by the legal establishment.

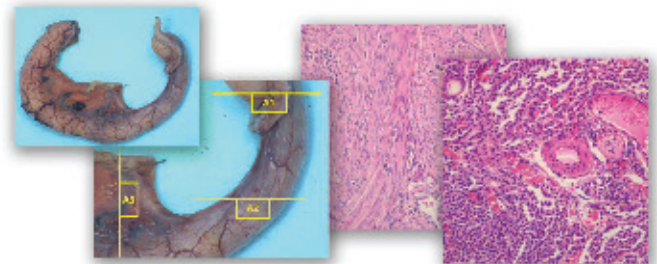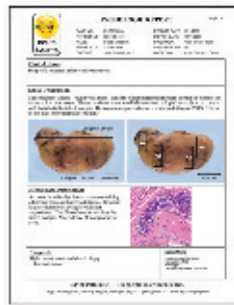

PathStand™ | The Smart Choice

Dimensions

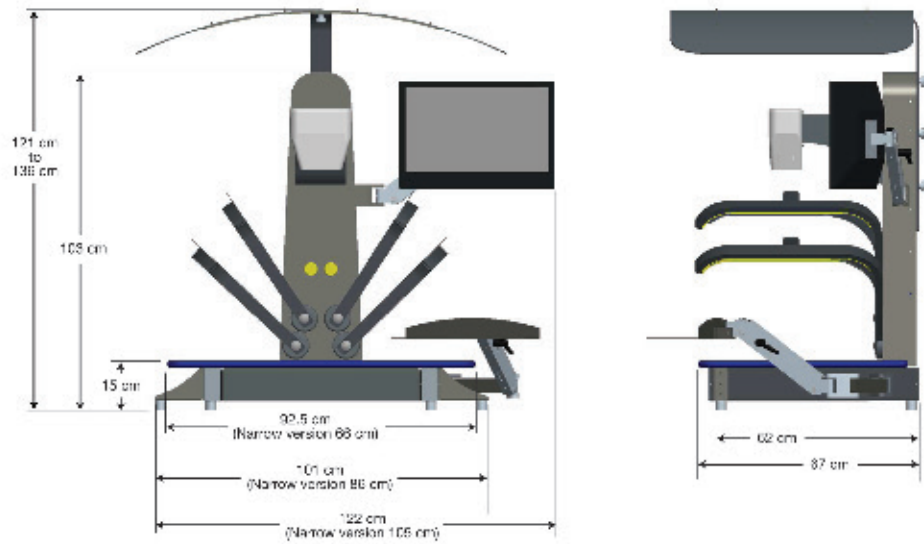

Easy to Implement

SPOT Imaging Solutions has developed turnkey solutions with configuration options to support your facility's unique hardware and software infrastructure. The result is an optimized work environment that is easily implemented into your existing infrastructure.



"The installation was quick and we were imaging the next day!"

*Satisfied Pathology Assistant,
Canada*

We Support You!

You are never alone with PathSuite InSpect™ technical service. Our 47 year tradition of standing behind our equipment ensures your uptime performance with on-line, on-phone or on-site assistance.



"We had a different system that was sitting in a box because we couldn't use it, PathSuite is so easy to use... Thanks SPOT Imaging!"

Medical Center, MI

Pathology Solutions to Support Your Whole Team

Microscope Imaging

- Routine Anatomical Pathology
- Frozen Section Consults
- Needle Biopsy Consults
- Tumor Board Presentations
- Teaching and Training

Macro Imaging

- Routine Gross Dissection
- Frozen Section Consults
- Teaching and Training

PathStand™ | Industry Leading Performance

Administrators...

"...the savings we achieve with the PathSuite software makes it a no brainer investment justification. We feel that we are now a state of the art pathology operation..."

University, NC

Users...

"So intuitive and easy to use! ... now that we have SPOT, we get our work done quickly and without frustration! It works great!"

University, GA

IT/Technical Support...

"...SPOT stands behind their products with on line support, repair services and reliable products ... They've been there for us!"

Medical Center, PA

Some of our Customer Locations...

University of Miami	Children's Hospital LA	Appalachian Reg Healthcare System
West Virginia University Hospital	Carillion Canada	Tricore Reference Laboratories
University of Georgia	Central Hospitals of Montreal	Temple University Hospital
Laval University-Eastern Quebec Telepathology	Cedars-Sinai Medical Center	Duke University
Broward Health	Maine Medical Center	Kaiser Hospital
California Pacific Medical Center	The University of Florida	Novant Health
Carolinas HealthCare System	Reading Hospital	Nycom Inc.
Orlando Regional Health Systems	Defense Finance & Acct. Service/DFAS Rome	Virginia Hospital Center
Texas Children's Hospital	Children's Healthcare of Atlanta	University of Pittsburgh Medical Center
Engineering Systems Inc.	Presbyterian Healthcare Services of NM	University of Pennsylvania Health System
Michigan State University	Robert Wood Johnson University Hospital	Alberta Provincial Healthcare

Call SPOT Imaging Solutions for experienced advice on how to dramatically improve your pathology lab with SPOT imaging and telepathology



6540 Burroughs Avenue • Sterling Heights • Michigan 48314-2133 USA
Toll Free USA: **1-866-604-SPOT** • Phone: **+1 586-731-6000** • Fax: **+1 586-731-6469**
Web: **www.spotimaging.com**